



positioning
patients for life.®

NEW RapidHeat™ Oven





RapidHeat™ Hot Air Oven for Aquaplast RT™ and Fibreplast™

The RapidHeat™ Oven quickly reaches working temperatures and provides an easy to use method to heat Aquaplast RT™ and Fibreplast™ thermoplastic masks as well as Aquaplast RT™ Custom Bolus™. The unique design allows the hot air to envelop the thermoplastic to soften quickly and evenly. This reduces the overall mask making time. The RapidHeat™ Oven features an easy to read digital timer and temperature control. The two sliding racks boast a non-stick coating which reduces the potential for sticking of the thermoplastic. The large interior of the oven has the ability to heat multiple masks simultaneously. For added safety the RapidHeat™ Oven automatically shuts off after 2 hours of inactivity.

Additional Features:

- Minimizes cross contamination
- Non-stick racks for easy cleaning
- Easy to read 7.5 cm digital display

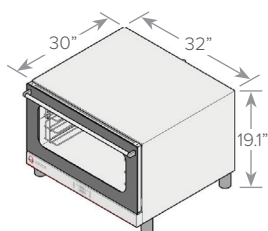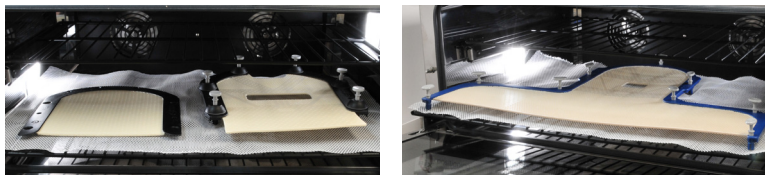

Outer dimensions: 32" w x 30" d x 19.1" h (inches)
81.3w x 76.2d x 48.5h cm

Inner dimensions: 27" w x 21" d x 10" h (inches)
68.6w x 53.3d x 25.4h cm


Features & Accessories



Digital Display

RapidHeat™ has a large 7.5 cm high digital display that features a Quick Start Guide to help easily navigate operating the oven. The digital temperature control can be used in a programmable or manual timer mode, increasing work flow efficiencies. The arrow indicators adjust time/temperature and the Start/Stop function allows you to easily preheat the oven while prepping your patient.

No more need for external timers, the digital display includes a timer that counts down from the designated time and alerts you when ready. The time can be set accordingly for different masks and Aquaplast RT™ Custom Bolus™.

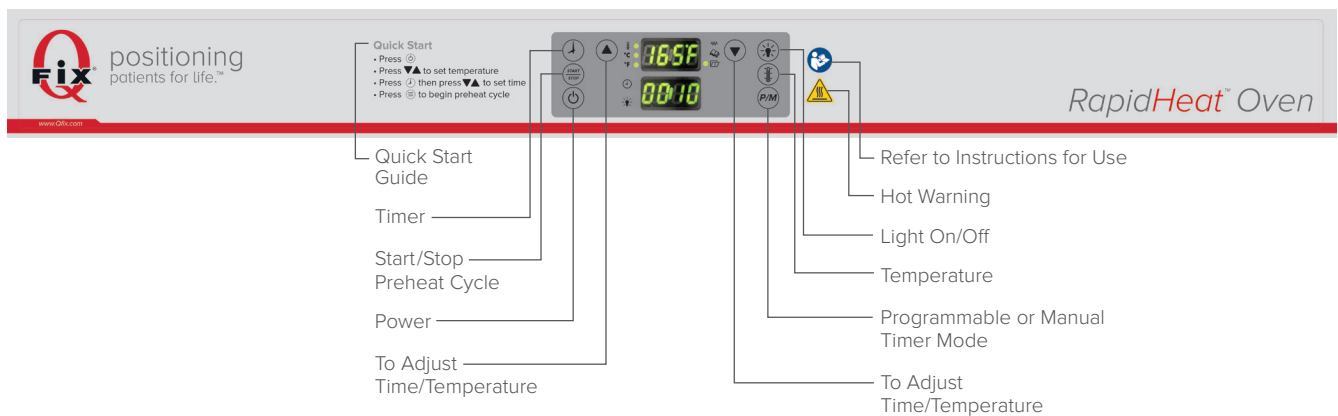

Thermoplastics Prep Cart

The easy to clean stainless steel Thermoplastic Prep Cart features two shelves for additional storage and wheels for easy transport. It is ideal for organizing, storing and transporting the RapidHeat™ Oven, Aquaplast RT™, Fibreplast™ and accessories.



Thermoplastics
Prep Cart
RT-2111B

RT-2075-110	RapidHeat™ Oven, 110 Volts
RT-2075-220	RapidHeat™ Oven, 220 Volts
RT-2111B	Thermoplastics Prep Cart



Qfix
Avondale, PA
+1 610.268.0585
sales@Qfix.com
www.Qfix.com