



# ONE Simple Solution for a Variety of Treatment Setups

## NEW kVue™ ONE Couch Top

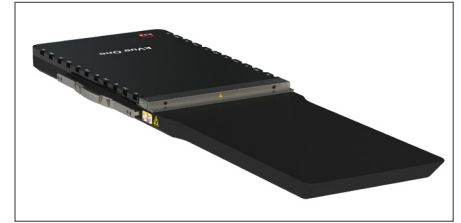

Without moving rails

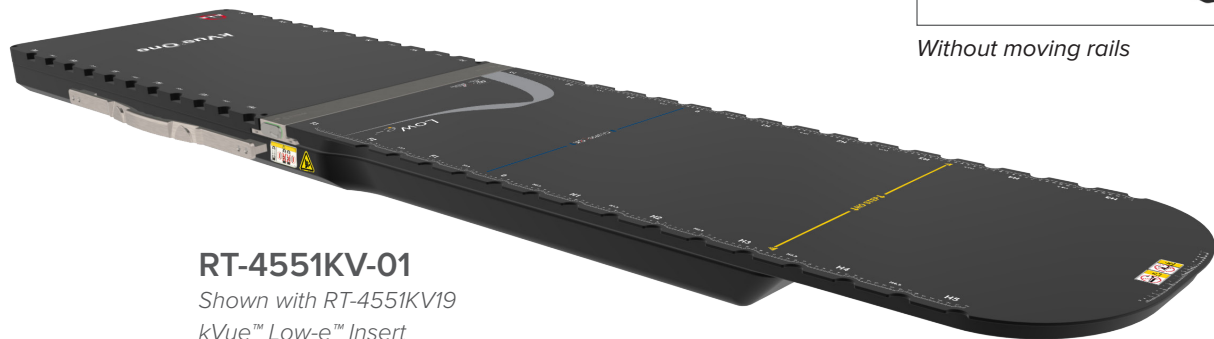

**RT-4551KV-01**


*Shown with RT-4551KV19  
kVue™ Low-e™ Insert*

The kVue™ One Couch Top is a one-piece modular design with a homogeneous composite base and no moving rails. It offers the ability to interchange treatment specific kVue™ inserts that mount directly into the top surface of the couch through the One Touch Latch™. The light weight kVue™ One is compatible with versions of 6 Degrees of Freedom (DOF) linear accelerators and robotic bases. (5 mm WET without kVue™ Insert)

## NEW kVue™ Low-e™ Insert


**Low**  
*e-*™

CALYPSO  OK

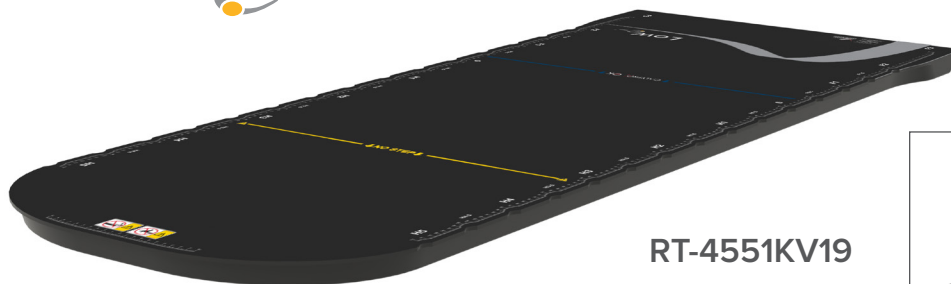


**RT-4551KV19**

*Shown with kVue™ SBRT Bridges  
RT-4558-B1 (small),  
RT-4558-B2 (large)*

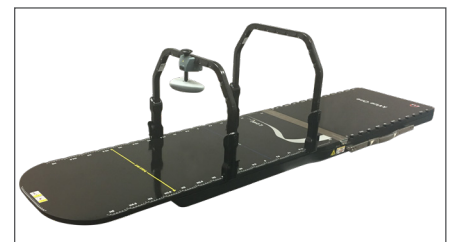

The kVue™ Low-e Insert is a low electron generating, homogeneous insert for enhanced imaging and treatment. It is the optimal choice for a wide range of radiotherapy treatments. (3.8 mm WET)

# SBRT Accessories

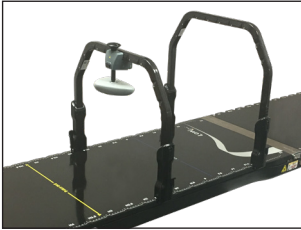

Shown with low profile compression paddle

## kVue™ SBRT Bridges

### RT-4558-B1 (small)

### RT-4558-B2 (large)

- Low profile bridges fit into most CT bore sizes
- Horizontal scale to aid in laser alignment
- Two sizes, Small (min: 30 cm & max: 38 cm)\* and Large (min: 39 cm & max: 48 cm)\* are completely adjustable in height
- \*Measurements are from the treatment surface to the underside of the bridge.



## SofTouch™ Knee Wedges (large & small) & Foot Block

### RT-4489KW (large)

### RT-4489KW-01 (small)

### & RT-4489FB

- Ergonomic design
- Indexes to locating bar in 35 mm increments
- Clip-on leg straps for added immobilization (included)
- Centering groove to aid in laser alignment



Shown with standard compression paddle

## Compression Paddle

### RT-4558KV4 (standard)

### RT-4558KV4-01 (low profile)


- Quick release button
- Tactile compression
- Incremental lateral placement along bridge



## ArmShuttle™ Elite

### RT-4546C

- Lightweight
- Ergonomic design
- MR Safe & CT compatible

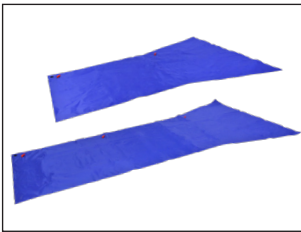


## 2-Chamber VacQfix™ Cushion

### RT-4517-0002

## 3-Chamber VacQfix™ Cushion

### RT-4517-0003


- Can be formed by one person
- Customizable full body immobilization
- Flared top for arm immobilization

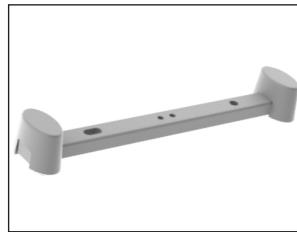

## VacQfix™ Indexer

### RT-4517-IND01

- Quick, easy vacuum cushion indexing
- Adaptable for a variety of setups

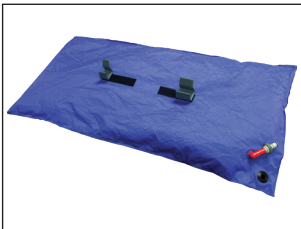

## Thigh Compression

## VacQfix™ Cushion

### RT-4517-4079VF15

- Complete leg immobilization
- Adjustable along bridge
- Snap-on attachment



## ZiFix™, Abdominal Thoracic

## Motion Control System

### RT-4558CB02

- Suppresses motion of lungs and liver for imaging & radiotherapy
- Compact, one-body design
- Three belt sizes

